
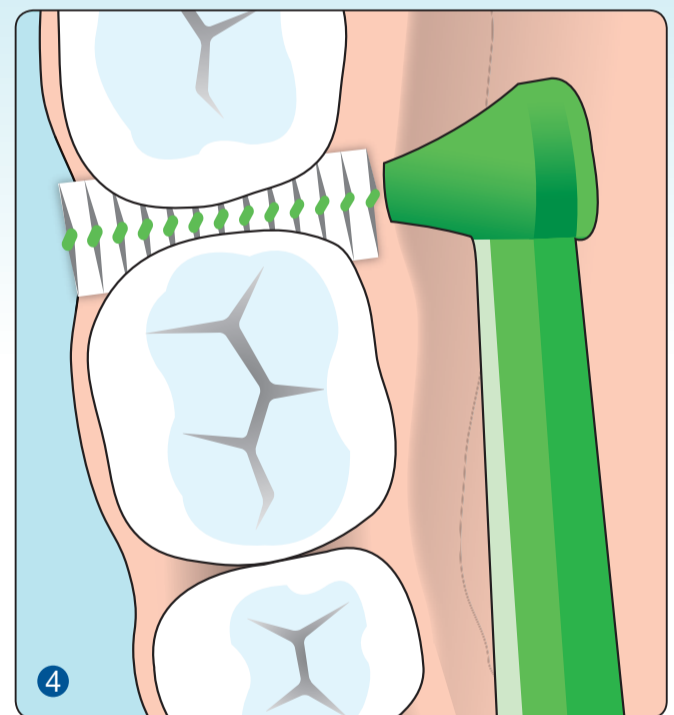
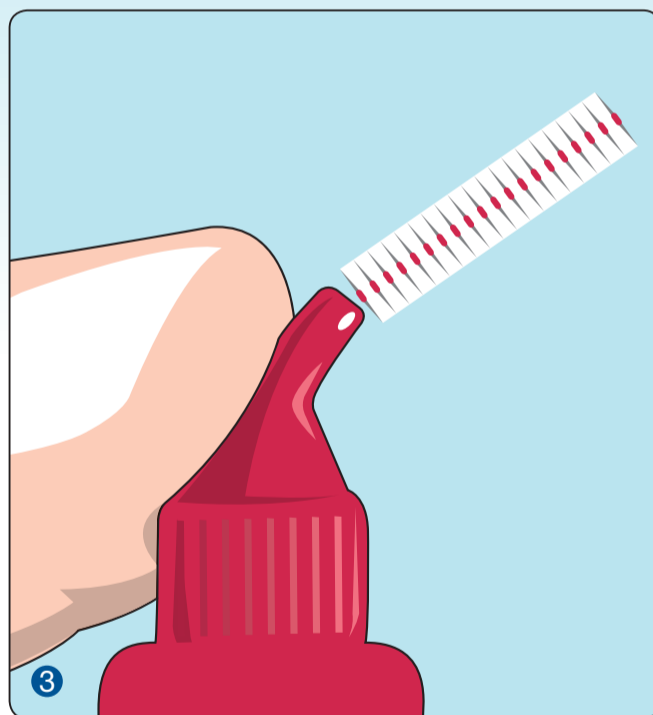
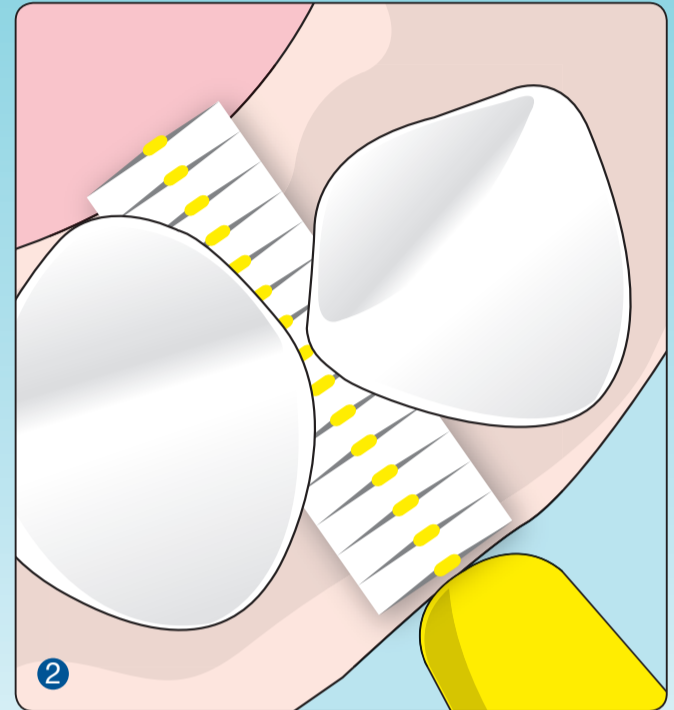


# Still only brushing 60% of your teeth?

TePe clean the 40% missed by toothbrushing alone.

- 1 Gently insert between the teeth at gum level, turning slightly. The correct size should fit easily but snugly.
- 2 Move brush back and forth a few times in each space. Never force, gently ease into a space.
- 3 Access may be improved by shaping the brush into a slight curve or by using our long handle brush – TePe Angle™.
- 4 TePe Angle™ is ideal for cleaning difficult-to-reach areas such as the back of the mouth.

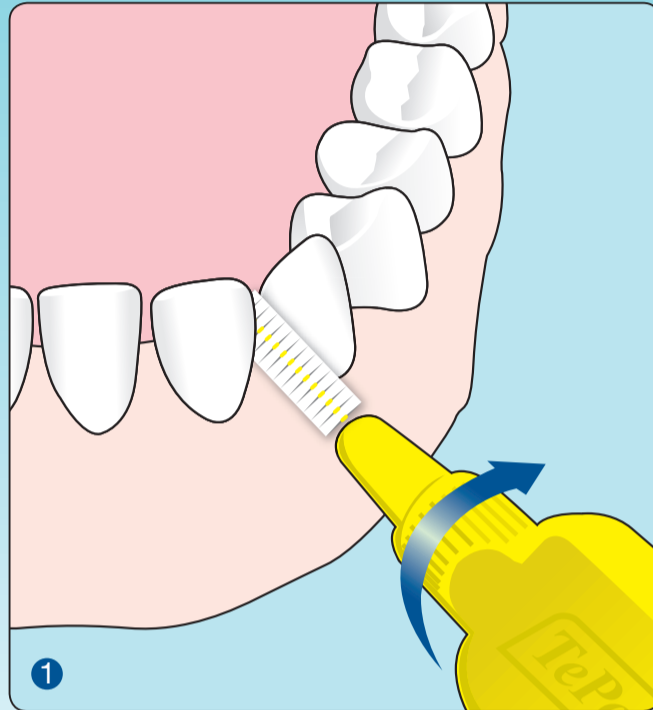

- ✓ Use daily
- ✓ Always rinse the brush after use
- ✓ Discard when filaments are worn

Ask your Dental Care Professional for help in selecting the correct brush size(s) for you.

TePe daily  
for best results!



TePe Original Interdental Brushes



TePe Angle™ Interdental Brushes

For reaching  
the back teeth



[www.tepe.co.uk](http://www.tepe.co.uk) TePeUK | [www.tepe.ie](http://www.tepe.ie)

